



*Tablets
Capsules
LCT's
Softgels*

Up to 70,000 PPH

Delta UV Laser

One or Two-Sided UV Laser
Single-Lane Carrier Link Feed System

Delta UV Laser

UV Laser Marking for Tablets & Capsules

The Delta UV Laser writer can laser write an alphanumeric, logo or trademark on a tablet, soft gelatin capsule, or caplet containing titanium dioxide (TiO2). The fully adjustable hopper allows the product to flow smoothly to the laser station and can laser mark 50,000 to 70,000 pieces per hour.

Features

- Single or two-sided laser marking
- Adjustable hopper for maximum fill rate
- Servo drive system for precision print registration
- Quick change tooling
- Swingarm 15" touchscreen HMI
- Wi-Fi connectivity
- Batch recording
- 21 CFR Part 11 Compliant
- Swing arm touch screen HMI
- CE and UL certified

Capabilities


- Output Up to 70,000 / Hour
- Single or Double Sided Marking with independent graphics for each side
- Mark alpha-numeric text, vector graphics, logos and bar codes
- Graphics editing software provided
- Adjustable parameters to control logo size, line width and etch depth
- Mark Bi-Color products with logo oriented on selected color side
- Orient logos based on product feature, such as tablet bisect line

Options

- Vision inspection system with single piece rejection
- Validation document packages available
- Floor level infeed bulk conveyor

Safety & Regulatory Compliance

- | | |
|-----------------------------|--------------------|
| • cGMP for Pharmaceuticals: | • NFPA 79 |
| FDA 21 CFR Compliance: | • IEC / EN 60204-1 |
| - Part 210 | • ISO 13849-1 |
| - Part 211 | • CRC, c 870 |
| - Part 11 | • CDRH Compliant |
| • CE Compliant | |

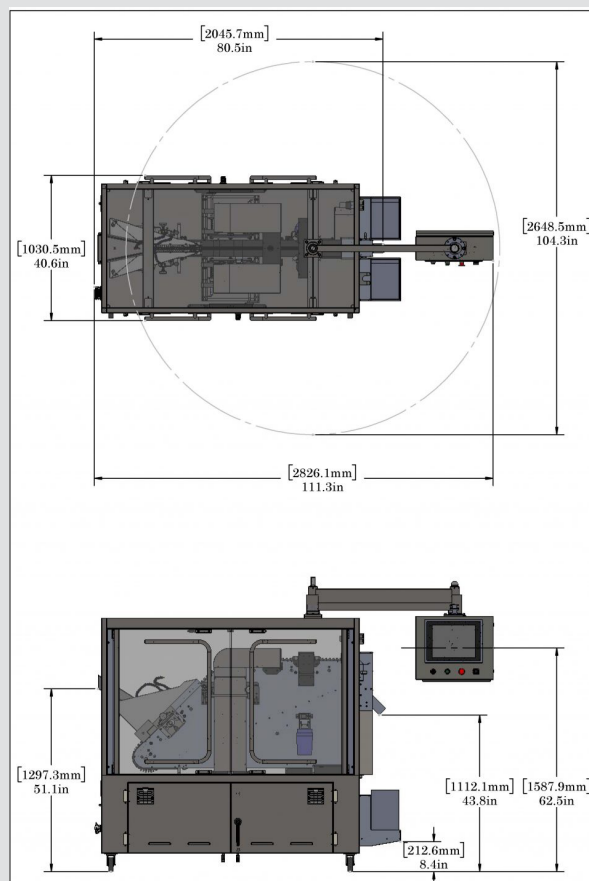

Specifications

H	2150mm / 84.6in
H1	1264.3mm / 49.8in
H2	1578.9mm / 62.2in
W1	829.2mm / 32.6in
W2	2416.4mm / 95.1in
L1	1663.5mm / 65.5in
L2	2328.5mm / 91.7in
POWER	240 VAC / 20 amp / 1 phase
AIR PSI	Air 90 psi @ 5 scfm



Made in the USA