



**Tablets
Capsules
LCT's
Softgels**

Up to 300,000 PPH

IBM C02 Laser

**One Sided CO2 Laser Writer
Bulk Hopper & Floor Feeder**

IBM C02 Laser

C02 Laser Marking for Tablets & Capsules

The IBM C02 Laser Writer is our high-speed, large output laser writing machine. With outputs over 300,000 pieces per hour, the IBM marks all soft gelatin capsules and caplets with ease. This flatbed machine comes with a bulk hopper and flood feeder for maximum fill rates and can print one side of a tablet.

Features

- Single sided laser marking
- Adjustable hopper for maximum fill rate
- Servo drive system for precision print registration
- Quick change tooling
- Swingarm 15" touchscreen HMI
- Wi-Fi connectivity
- Batch recording
- 21 CFR Part 11 Compliant
- CE and UL certified

Capabilities


- Output Up to 300,000 / Hour
- Single Sided Marking
- Mark alpha-numeric text, vector graphics, logos and bar codes
- Graphics editing software provided
- Adjustable parameters to control logo size, line width and etch depth
- Mark Bi-Color products with logo oriented on selected color side
- Orient logos based on product feature, such as tablet bisect line

Options

- Vision inspection system with single piece rejection
- Validation document packages available
- Floor level infeed bulk conveyor

Safety & Regulatory Compliance

- | | |
|-----------------------------|--------------------|
| • cGMP for Pharmaceuticals: | • NFPA 79 |
| FDA 21 CFR Compliance: | • IEC / EN 60204-1 |
| - Part 210 | • ISO 13849-1 |
| - Part 211 | • CRC, c 870 |
| - Part 11 | • CDRH Compliant |
| • CE Compliant | |

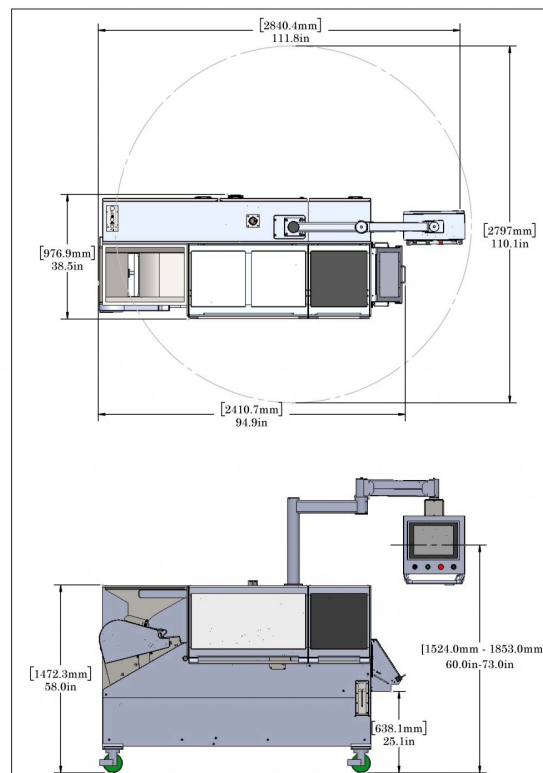

Specifications



Made in the USA

H	2314.4mm / 91.1in
H1	1480.3mm / 58.3in
H2	1524mm - 1853mm / 60.0in - 73.0in
W1	955.5mm / 37.6in
W2	2740.6mm / 107.9in
L1	2832.9mm / 111.5in
L2	2352.8mm / 92.6in
POWER	240 VAC / 30 amp / single phase
AIR PSI	Air 90 psi @ 5 scfm